

# Urilyzer® 500 Pro



The Urilyzer® 500 Pro is a semi-automated urine test strip analyzer with enhanced functionality to the proven CombiScan® 500.

The high throughput of up to 500 tests/hour in combination with state of the art operator and data management functions provides the ideal system for mid-sized to bigger laboratories and for hospitals.

## Advanced Solution in Urine Diagnostics



The Urilyzer® 500 Pro provides enhanced operator management and connectivity capabilities to meet the growing demands for compliance management and data capture.

In combination with the urine test strip CombiScreen® 11SYS and 11SYS PLUS the system offers uncompromising ease of use, high quality and efficiency in a compact format.

- Easy to use
- Data Management
- Quality Control Management
- Operator Management

### Easy to use

- Automatic strip detection (autostart)
- Easy to use software
- Color touchscreen operation
- High throughput of up to 500 tests/hour
- Microscopy flag (sediment recommendation)
- Standby function
- Minimal maintenance and convenient cleaning
- Continuous loading system

### Data Management


- Flexible, customized testing and reporting options
- Advanced data entry (e.g. color, clarity, lot, expiry)
- Worklist transmission through LIS2 and HL7
- Optional barcode reader and external keyboard

### Quality Control Management

- Automated test strip check
- Automated QC analysis (QC Reminder, QC lockout function)
- Lot number and expiration date check

### Operator Management

- Multiple authorization levels
- Programmable operator management provides customized access levels & prevents unauthorized use
- Settings and operator list import/export via USB stick

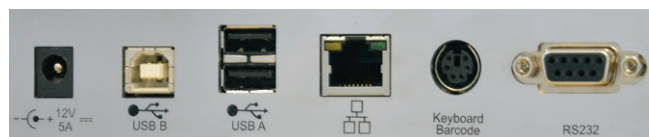

## Specifications

Type	Semi-automated urine test strip analyzer with continuous loading
Measurement technology	Reflectance photometer with 4 wavelengths: 505, 530, 620, 660 nm
Parameters	Bilirubin, Urobilinogen, Ketones, Ascorbic Acid, Glucose, Protein (Albumin), Blood (Hemoglobin), pH, Nitrite, Leucocytes, Specific Gravity
Throughput	Up to 500 tests/hour
Data storage	Patient database: 5000 tests / QC database: 5000 tests
Display	5.7" LCD VGA capacitive touchscreen
Interfaces	Serial RS232, USB Type A, USB Type B, PS2 (external keyboard, barcode reader), Ethernet
Dimensions	300 x 300 x 180 mm (WxDxH)
Weight	6.5 kg
Power supply	12 V DC / 5 A
Operating environment	Temperature: +15°C to +32°C Relative humidity (non-condensing): 20% to 85%
Printer	Built-in thermal printer
Barcode Reader	External; Barcode types: Code 39, Code 128, Interleaved 2/5, EAN-8
Protocols	LIS2 (ASTM+), HL7

Urine strip analyzer Urilyzer® 500 Pro, Art.-No.: UL0500Pro

Distributor information



### Analyticon Biotechnologies AG

Am Muehlenberg 10  
35104 Lichtenfels - Germany  
Phone: +49 6454 7991-0  
info@analyticon-diagnostics.com  
www.analyticon-diagnostics.com